FreeStyle *LibreLink*

**User's Manual
Manuel d'utilisation**

English

Français

X

App Symbols**Important Information**

Indications for Use

FreeStyle LibreLink Overview

Home Screen

FreeStyle Libre Sensor Kit

App Setup**Applying Your Sensor****Starting Your Sensor****Checking Your Glucose****Understanding Your Glucose Readings****Adding Notes****Reviewing Your History**

Logbook

Other History Options

Removing Your Sensor**Replacing Your Sensor****Setting Reminders****Settings and Other Options
in the Main Menu****Activities****Maintenance and Disposal****Troubleshooting**

Problems at the Sensor Application Site

Problems Starting Your Sensor
or Receiving Sensor Readings**Customer Service****Labelling Symbols and Definitions****Electromagnetic Compatibility****Performance Characteristics**

User's Manual

App Symbols

Sensor may be inaccurate.
Check blood glucose with
a test strip before making
any treatment decisions.

Direction your glucose is
going. See [Understanding
Your Glucose Readings](#) for
more information.

Scan button

Caution

Add/edit notes

Food note

Insulin (Rapid or Long-
acting) note

Exercise note

Time change

Sensor too cold

Sensor too hot

App icon

Multiple/Custom notes

Share report

Additional information

Main menu

Calendar

Important Information

Indications for Use

FreeStyle LibreLink app ("App") when used with a FreeStyle Libre Flash Glucose Monitoring System Sensor ("Sensor") is indicated for measuring interstitial fluid glucose levels in adults aged 18 years and older with diabetes mellitus. The App and Sensor are designed to replace blood glucose testing in the self-management of diabetes, including dosing of insulin. Treatment decisions should not be based on real-time Sensor glucose readings alone and instead should consider all the information on the results screen.

WARNING: If you are using FreeStyle LibreLink, you must also have access to a blood glucose monitoring system as the App does not provide one.

CAUTION: FreeStyle LibreLink installed on a smartphone is intended for use by a single person. It must not be used by more than one person due to the risk of misinterpreting glucose information.

Additional Safety Information

FreeStyle LibreLink and FreeStyle Libre Flash Glucose Monitoring System Readers ("Readers") do not share data. For complete information on a device, be sure to scan your Sensor every 8 hours with that device; otherwise, your reports will not include all your data.

Security Information

- You are responsible for properly securing and managing your smartphone. If you suspect an adverse cybersecurity event related to FreeStyle LibreLink, contact Customer Service.
- FreeStyle LibreLink is not intended for use on a smartphone that has been altered or customised to remove, replace, or circumvent the manufacturer's approved configuration or use restriction, or that otherwise violates the manufacturer's warranty.

The following Contraindication, Warnings, and other safety information apply to the Sensor, when used with FreeStyle LibreLink.

CONTRAINDICATION

MRI/X-ray/CT: Remove the Sensor you are wearing if you have a medical appointment that includes strong magnetic or electromagnetic radiation e.g. an X-ray, MRI (Magnetic Resonance Imaging), or CT (Computed Tomography) scan. Apply a new Sensor after the appointment.

WARNING:

- **Choking Hazard:** The Sensor contains small parts that may be dangerous if swallowed.
- When your glucose is changing quickly (more than 0.1 mmol/L per minute), interstitial fluid glucose levels measured by the Sensor may not accurately match blood glucose levels. When your glucose is rising quickly or falling quickly, check the Sensor glucose readings by doing a blood glucose test on your finger with a blood glucose meter.
- Do a blood glucose test on your finger with a blood glucose meter if you get a Low Glucose or Glucose Going Low message.
- Do a blood glucose test on your finger with a blood glucose meter if you do not see a Current Glucose number or Glucose Trend Arrow.
- Do not ignore symptoms that may be due to low or high blood glucose. If you have symptoms that do not match the Sensor glucose reading or think that your reading may be inaccurate, check the reading by doing a blood glucose test on your finger with a blood glucose meter. If you have symptoms that do not match your glucose readings, consult your health care professional. **If you have been diagnosed with hypoglycaemic unawareness and do not experience symptoms of hypoglycaemia, ask your health care professional when you should perform a blood glucose test to confirm your glucose. The System**

has not been evaluated for use in patients with hypoglycaemic unawareness and will not automatically alert you of a hypoglycaemic event without you scanning your Sensor. For example, the System does not have an alarm that can alert or wake you when you are sleeping in the case of low glucose.

- When you see the  symbol, you must check your blood glucose with a blood glucose meter before making any treatment decisions. Sensor readings may not accurately reflect blood glucose levels.

CAUTION:

- On rare occasions, you may get inaccurate Sensor glucose readings. If you believe your readings are not correct or are inconsistent with how you feel, perform a blood glucose test on your finger to confirm your glucose and check to make sure your Sensor has not come loose. If the problem continues or if your Sensor is coming loose, remove the current Sensor and apply a new one.
- Intense exercise may cause your Sensor to loosen due to sweat or movement of the Sensor. If your Sensor comes loose, you may get no readings or unreliable readings, which may not match how you feel. Follow the instructions to select an appropriate application site.
- The System uses all available glucose data to give you readings so you should scan your Sensor at least once every 8 hours for the most accurate performance. Scanning less frequently may result in decreased performance. If you are using both the App and Reader with the same Sensor, be sure to scan frequently with both devices.
- Some individuals may be sensitive to the adhesive that keeps the Sensor attached to the skin. If you notice significant skin irritation around or under your Sensor, remove the Sensor and stop using the Sensor. Do not try to treat any skin irritation with lotions/ creams/ ointments/ sprays/ barrier patches in order to continue wearing the Sensor. Contact your health care professional before continuing to use the Sensor.
- When using the System to replace blood glucose testing for making diabetes treatment decisions, including insulin dosing, you should have a good understanding of how to self-manage your diabetes, as determined by your health care professional.
- Performance of the Sensor when used with other implanted medical devices, such as pacemakers, has not been evaluated.

- Do not reuse Sensors. The Sensor has been designed so that it cannot be reused. Not suitable for re-sterilisation.
- Severe dehydration and excessive water loss may cause inaccurate Sensor results. If you think you are dehydrated, consult your health care professional immediately.
- The Sensor Pack and Sensor Applicator are packaged as a set and have the same Sensor code. Check that the Sensor codes match before using your Sensor Pack and Sensor Applicator. Sensor Packs and Sensor Applicators with the same Sensor code should be used together or your Sensor glucose readings may be incorrect.

Additional Safety Information

- Glucose levels in the interstitial fluid can be different from blood glucose levels and may mean that Sensor glucose readings are different from blood glucose. You may notice this difference during times when your blood glucose is changing quickly, for example after eating, taking insulin, or exercising.
- Glucose Going Low and Glucose Going High messages may mean that your actual blood glucose is already <3.9 mmol/L or >13.3 mmol/L. Consider this before making an activity or treatment decision.
- If you are going to participate in activities that could result in harm to yourself or others in the event of a severe hypoglycaemic episode (e.g. driving a motor vehicle without following the Canadian Diabetes Association Recommendations for Private and Commercial Drivers) do not rely on Glucose Messages (High Glucose, Low Glucose, Glucose Going High and Glucose Going Low) alone. Talk to your healthcare professional about checking your glucose before and while driving.
- **How to perform a Sensor accuracy check:** Check your Sensor is working properly by using your blood glucose meter. To check your Sensor, carry out the following steps:
 1. Scan your Sensor. Make sure your glucose is changing slowly →.
 2. Perform a blood glucose test. Your blood glucose result must be above 4.4 mmol/L in order to perform the check.
 3. Compare the blood glucose result to the Sensor reading. The Sensor is performing correctly when the blood glucose result is within 20% of the Sensor reading. For example, when the Sensor reads 10.0 mmol/L, the blood glucose

result would be expected to be between 8.0 and 12.0 mmol/L.

- **Interfering Substances:** Taking ascorbic acid (vitamin C) while wearing the Sensor may falsely raise your Sensor glucose readings. Taking salicylic acid (used in some pain relievers such as aspirin and some skin care products) may slightly lower your Sensor glucose readings. The level of inaccuracy depends on the amount of the interfering substance active in your body.
- Store the Sensor Kit between 4°C and 25°C. While you don't need to keep your Sensor Kit in a refrigerator, you can as long as the refrigerator is between 4°C and 25°C.
- The Sensor has not been evaluated for use in pregnant women, persons on dialysis, or people less than 18 years of age.
- The Sensor Pack is sterile unless opened or damaged.
- Your Sensor has been tested to withstand immersion into one metre (3 ft) of water for up to 30 minutes.
- Do not freeze the Sensor. Do not use if expiry date has passed.

FreeStyle LibreLink Overview

IMPORTANT: Read all of the information in this User's Manual before using FreeStyle LibreLink with a Sensor. Refer to your iPhone instructions for use for how to use your iPhone. If you are using a Reader, refer to the User's Manual in the Reader Kit.


FreeStyle LibreLink is available for download from the App Store. When you're ready to start using FreeStyle LibreLink, you'll prepare and apply a Sensor to the back of your upper arm. You can then use the App to get glucose readings from the Sensor and store your glucose history and notes. The Sensor can be worn on your body for up to 14 days.

Note:

- The Sensor comes in the FreeStyle Libre Sensor kit. See [FreeStyle Libre Sensor Kit](#).
- Go to www.FreeStyleLibre.com for smartphone requirements and compatibility. Please keep in mind that the ease of scanning a Sensor may vary between devices.

Home Screen

The Home Screen gives you access to information about your glucose and the App. To return to the Home Screen from another screen, go to the Main Menu and tap **Home**.

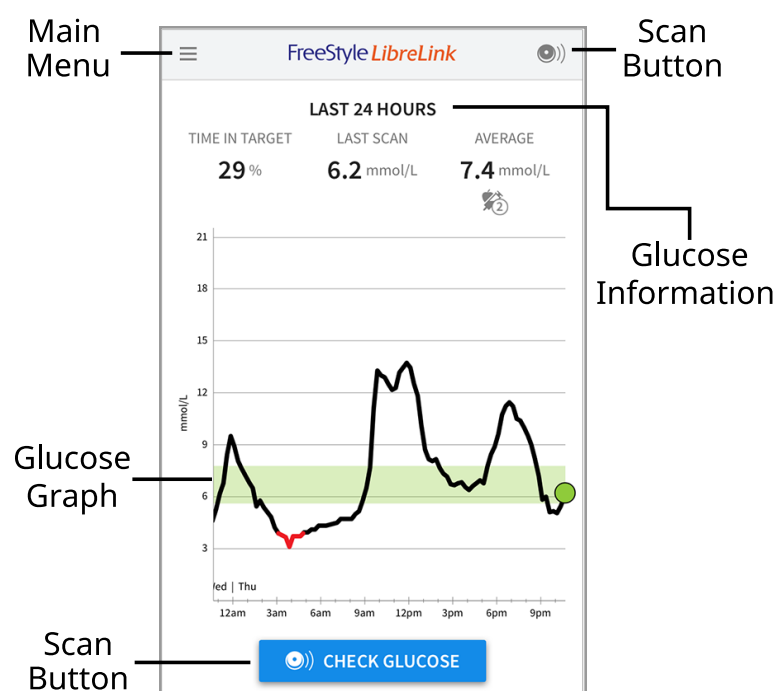Main Menu - Tap to access the Home Screen, Logbook, other history options, and the Share option. You can also access Settings, Help, and other information.

Glucose Graph - Graph of your stored Sensor glucose readings.

Scan Button - Tap when you are ready to scan your Sensor. You can either tap the blue box on the Home Screen or  at the top right.

Glucose Information - Your Time In Target, information about your last scan, and average glucose for the last 24 hours.

FreeStyle Libre Sensor Kit

The FreeStyle Libre Sensor Kit includes:


- Sensor Pack
- Sensor Applicator
- Alcohol wipe
- Product insert

When opening your kit, check that the contents are undamaged and that you have all parts listed. If any parts are missing or damaged, contact Customer Service. The Sensor (only visible after applied) is initially in two parts: one part is in the Sensor Pack and the other part is in the Sensor Applicator. Once prepared and applied to your body, the Sensor measures your glucose using a small, flexible tip that inserts just under the skin.

Sensor Pack. Used with the Sensor Applicator to prepare the Sensor for use.

Sensor Applicator. Applies the Sensor to your body.

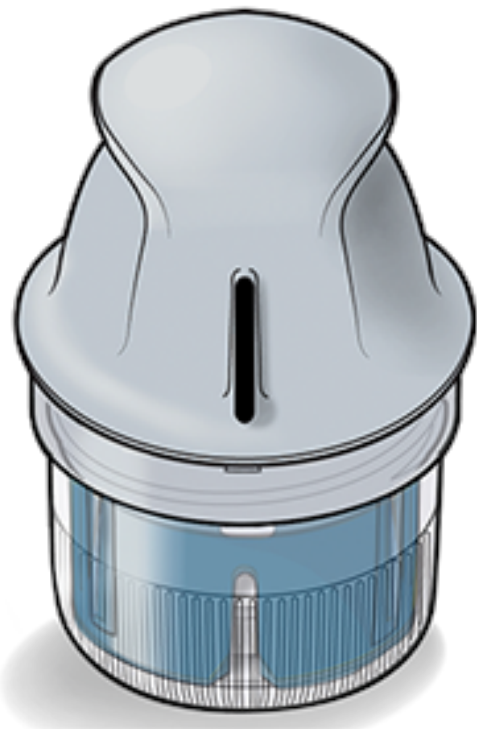App Setup

Before using the App for the first time, you must complete the setup.

1. Check that your iPhone is connected to a network (WiFi or mobile). You can then install FreeStyle LibreLink from the App store. Tap the App icon to open the App.

Note: You only need to be connected to a network for setup, using LibreView, and sharing with other apps. You do not need to be connected to scan a Sensor, add notes, or review your history in the App.

2. Swipe left to view some helpful tips or tap **GET STARTED NOW** at any point.
3. Confirm your country and tap **NEXT**.
4. You need a LibreView account to use the App. Follow onscreen instructions to review legal information and create a new account or login to your existing account.

LibreView Data Management Software is developed and distributed by Newyu, Inc. Use of FreeStyle LibreLink requires registration with LibreView, a service provided by Abbott and Newyu, Inc.

5. Confirm your glucose unit of measure and tap **NEXT**.
6. Set your Target Glucose Range and tap **NEXT**. Work with your health care professional to determine your Target Glucose Range. Your Target Glucose Range is displayed on glucose graphs in the App and is used to calculate your Time In Target.


7. Select how you count carbohydrates (in grams or portions) and tap **NEXT**. The carbohydrate unit will be used in any food notes you enter in the App.
8. The App now displays useful information about the My Glucose screen. Tap **NEXT** to view information about the Glucose Trend Arrow. Tap **NEXT** again.
9. Apply a new Sensor and then tap **NEXT**. Go to [Starting Your Sensor](#).

Note: If you need help applying your Sensor, tap **HOW TO APPLY A SENSOR** or go to [Applying Your Sensor](#).

Applying Your Sensor


CAUTION:

- The Sensor Pack and Sensor Applicator are packaged together as a set and have the same Sensor code. Check that the Sensor codes match before using your Sensor Pack and Sensor Applicator. Sensor Packs and Sensor Applicators with the same Sensor code should be used together or your Sensor glucose readings may be incorrect.

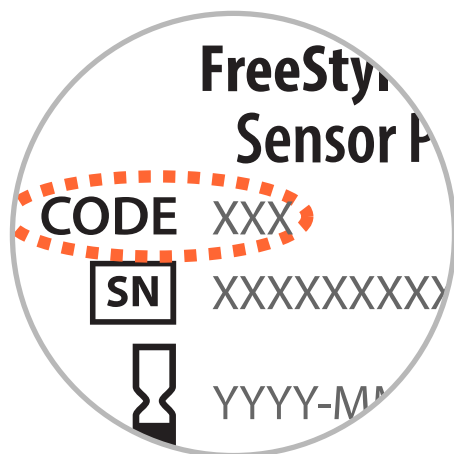- Intense exercise may cause your Sensor to loosen due to sweat or movement of the Sensor. If your Sensor comes loose, you may get no readings or unreliable readings, which may not match how you feel. Follow the instructions to select an appropriate application site.

-
1. Apply Sensors only on the back of your upper arm. Avoid areas with scars, moles, stretch marks or lumps. Select an area of skin that generally stays flat during your normal daily activities (no bending or folding). Choose a site that is at least 2.5 cm (1 inch) away from an insulin injection site. To prevent discomfort or skin irritation, you should select a different site other than the one most recently used.

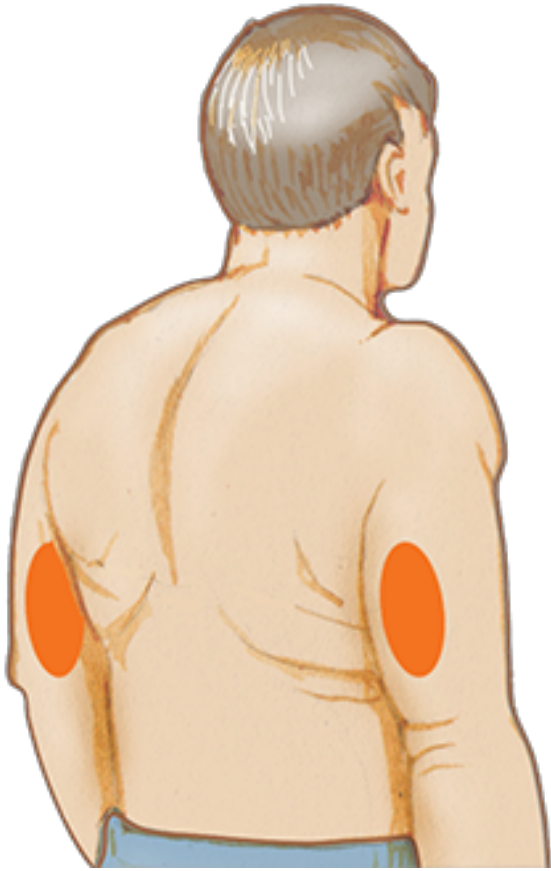2. Clean application site with an alcohol wipe and allow site to dry before proceeding. This helps the Sensor stay attached to your body.

Note: The area **MUST** be clean and dry, or the Sensor may not stick to the site.


3. Open the Sensor Pack by peeling the lid off completely. Unscrew the cap from the Sensor Applicator and set the cap aside.

CAUTION: Do NOT use if the Sensor Pack or Sensor Applicator seem to be damaged or already opened. Do NOT use if past expiry date.

4. Line up the dark mark on the Sensor Applicator with the dark mark on the Sensor Pack. On a hard surface, press down firmly on the Sensor Applicator until it comes to a stop.

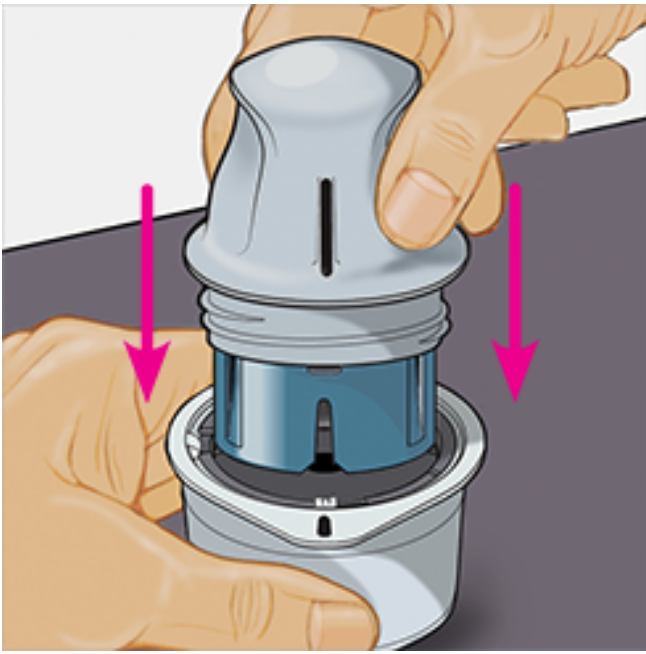
5. Lift the Sensor Applicator out of the Sensor Pack.


6. The Sensor Applicator is prepared and ready to apply the Sensor.

CAUTION: The Sensor Applicator now contains a needle. Do NOT touch inside the Sensor Applicator or put it back into the Sensor Pack.

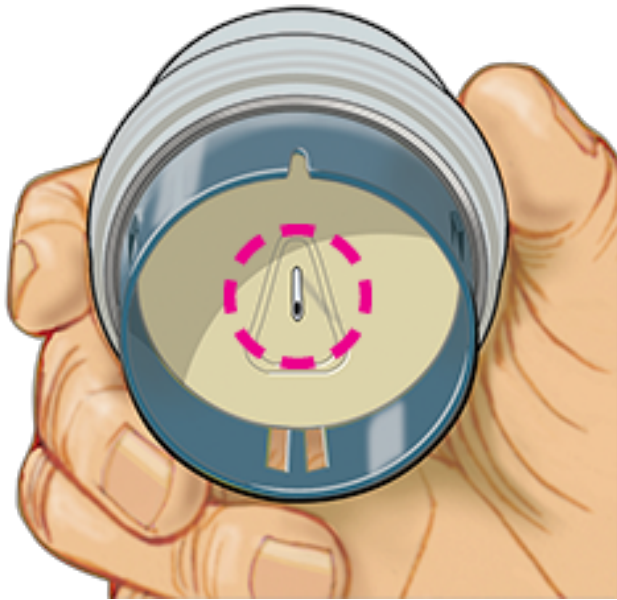
7. Place the Sensor Applicator over the prepared site and push down firmly to apply the Sensor to your body.

CAUTION: Do NOT push down on Sensor Applicator until placed over prepared site to prevent unintended results or injury.

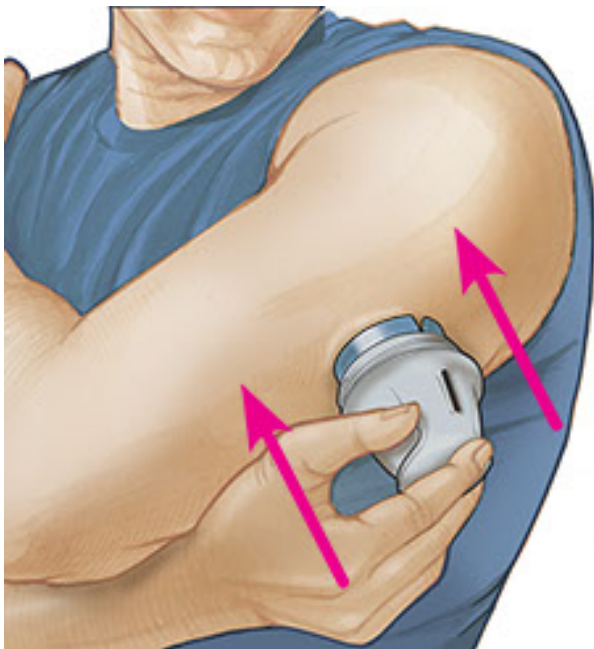8. Gently pull the Sensor Applicator away from your body. The Sensor should now be attached to your skin.

Note: Applying the Sensor may cause bruising or bleeding. If there is bleeding that does not stop, remove the Sensor, and apply a new one at a different site.

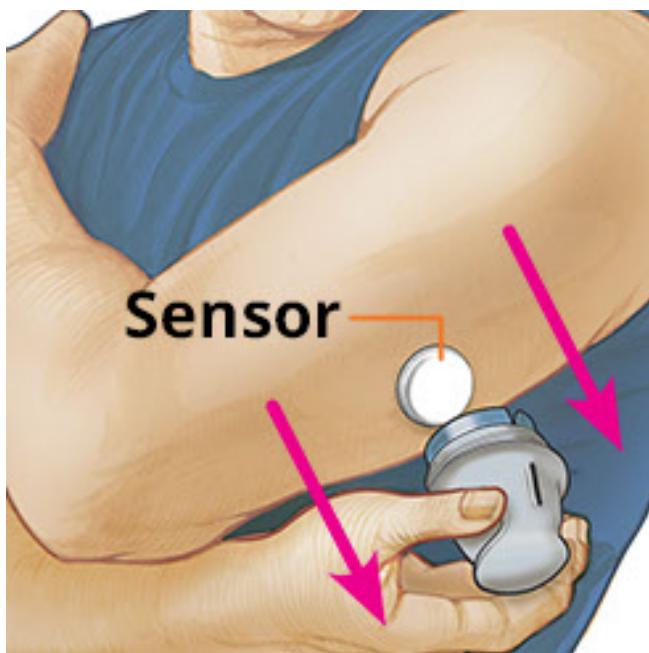9. Make sure Sensor is secure after application. Put the cap back on the Sensor Applicator. Discard the used Sensor Applicator and Sensor Pack according to local regulations.

Note: Tap **Help** in the Main Menu to access an in-app tutorial on applying a Sensor.

Starting Your Sensor

IMPORTANT:

- The App requires that your iPhone has date and time enabled to set automatically. You can check this in your iPhone settings.
- When using the App, you should keep your iPhone well charged and be sure you

have access to a blood glucose meter.

- When you scan your Sensor, you will receive a tone and vibration. If your iPhone's volume is turned off, you will not hear the tone.
 - The NFC (Near Field Communication) antenna is on the top edge of iPhone. Hold this area near your Sensor when you are scanning. You may need to adjust your scan distance based on what clothing you are wearing. In addition to proximity and orientation, other factors can affect NFC performance. For example, a bulky or metallic case can interfere with the NFC signal. Keep in mind that the ease of scanning a Sensor may vary between devices.
-

1. Tap the scan button .

Note:

- You can either tap the blue box on the Home Screen or  at the top right.
- If the Ready to Scan dialogue disappears, tap the scan button again.

NFC is now activated and your iPhone is ready to scan the Sensor.

2. Hold the top of your iPhone near the Sensor (this can be done over clothing). Do not move your iPhone until you hear a tone and/or feel a vibration. This completes the scan.

Note:

- If you need help, tap **HOW TO SCAN A SENSOR** to view an in-app tutorial. You can also access this later by going to the Main Menu and then tapping **Help**.
- If your Sensor is not successfully scanned, you may receive this Scan Error: "Your scan was unsuccessful. Tap the scan button and scan again."

See [Troubleshooting](#) for additional error messages.

3. The Sensor can be used to check your glucose after 60 minutes. While the Sensor is starting up you can navigate away from the App. If notifications are enabled, you will see a notification when the Sensor is ready.

Note:

- If you like, you can use a Sensor with both the App and the Reader. To do this, you must start the Sensor with the Reader first and then scan with the

App.


- Remember that FreeStyle LibreLink and Readers do not share data. For complete information on a device, be sure to scan your Sensor every 8 hours with that device; otherwise, your reports will not include all your data.

Checking Your Glucose

1. Open the App and tap the scan button .

Note: If the Ready to Scan dialogue disappears, tap the scan button  again.

2. Hold the top of your iPhone near the Sensor until you hear a tone and/or feel a vibration.
3. The My Glucose screen now displays your glucose reading. It includes your Current Glucose, a Glucose Trend Arrow indicating which way your glucose is going, and a graph of your current and stored glucose readings.

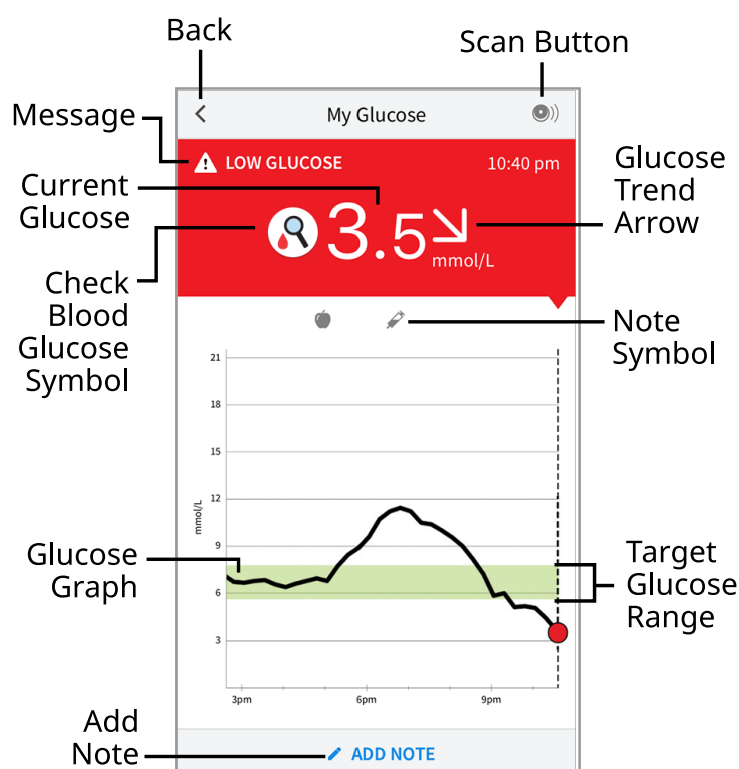Scan Button - Tap when you are ready to scan your Sensor.

Message - Tap for more information.

Check Blood Glucose Symbol - When you see this symbol , do a blood glucose test before making treatment decisions.

Back - Tap to return to Home screen.

Current Glucose - Glucose value from your latest scan.

Add Note - Tap to add notes to the glucose reading.

Glucose Trend Arrow - Direction your glucose is going.

Note Symbol - Tap to review notes you've entered.

Glucose Graph - Graph of your current and stored glucose readings.

Target Glucose Range - The graph shows your target glucose range.

Note:

- A Sensor can store up to 8 hours of glucose data, so scan it at least once every 8 hours to capture all of your available glucose data.
- The graph will scale to 27.8 mmol/L to accommodate glucose readings above 21 mmol/L.
- The  symbol may appear, indicating the smartphone's time was changed. Gaps in the graph may result or glucose readings may be hidden.
- Your current glucose value determines the background colour on the My Glucose screen:

Orange - High glucose (above 13.3 mmol/L)

Yellow - Between the Target Glucose Range and high or low glucose level

Green - Within the Target Glucose Range

Red - Low glucose (below 3.9 mmol/L)

Understanding Your Glucose Readings

Glucose Trend Arrow


The Glucose Trend Arrow gives you an indication of the direction your glucose is going.

Glucose is rising quickly
(more than 0.1 mmol/L
per minute)

Glucose is rising
(between 0.06 and
0.1 mmol/L per minute)

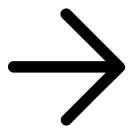Glucose is changing
slowly (less than
0.06 mmol/L per minute)

Glucose is falling
(between 0.06 and
0.1 mmol/L per minute)

Glucose is falling quickly
(more than 0.1 mmol/L
per minute)


Note: The Glucose Trend Arrow may not always appear with your reading.

Messages


Below are messages you may see with your glucose readings.

If your glucose reading is less than 3.9 mmol/L, projected to be less than 3.9 mmol/L, rapidly changing, or there is no number or trend arrow, you will see this symbol . You can touch the symbol for more information. Check your blood glucose on your finger with a test strip before making treatment decisions.


LO | HI: If **LO** appears, your reading is lower than 2.2 mmol/L. If **HI** appears, your reading is higher than 27.8 mmol/L. You can touch the  symbol for more information. Check your blood glucose on your finger with a test strip. If you get a second **LO** or **HI** result, contact your health care professional **immediately**.

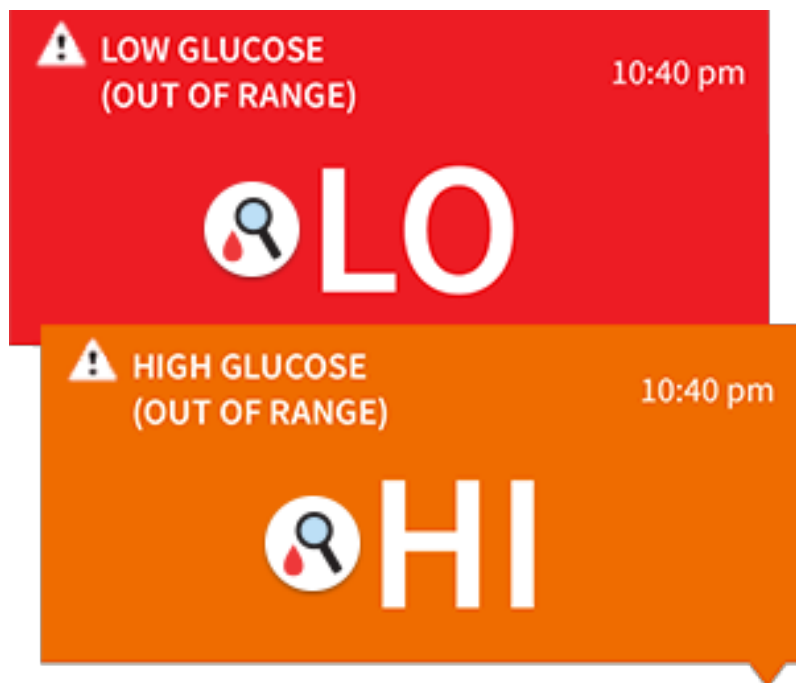Low Glucose | High Glucose: If your glucose is higher than 13.3 mmol/L or lower than 3.9 mmol/L, you will see a message on the screen. You can touch the  symbol for more information and set a reminder to check your glucose.

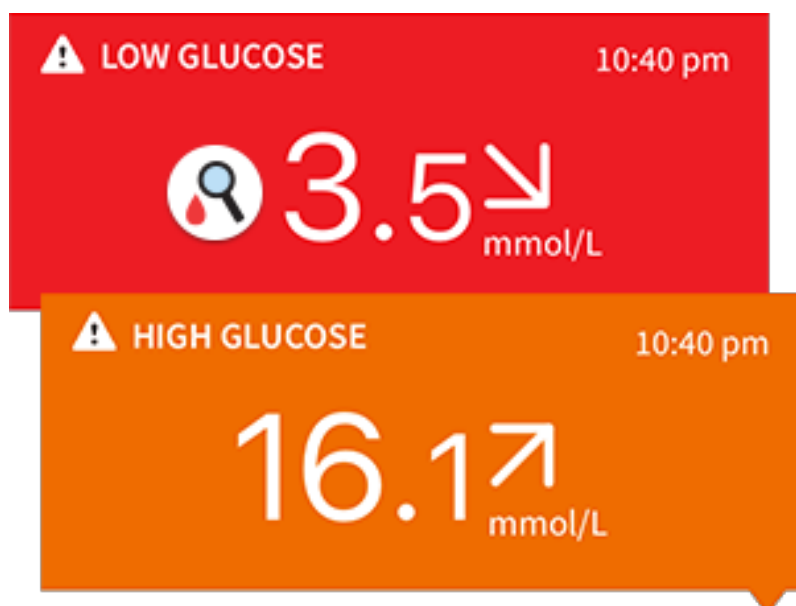Glucose Going Low | Glucose Going High: If your glucose is projected to be higher than 13.3 mmol/L or lower than 3.9 mmol/L within 15 minutes, you will see a message on the screen. The background colour corresponds to your current glucose value. You can touch the  symbol for more information and set a reminder to check your glucose.

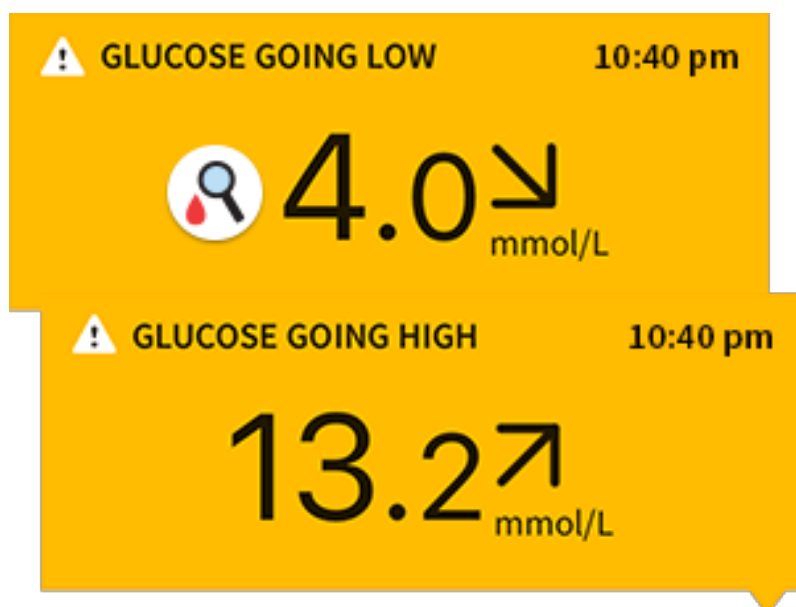Note: If you are not sure about a message or reading, contact your health care professional for information.

Adding Notes

Notes can be saved with your glucose readings to help you track food, insulin, and exercise. You can also add your own comment.

1. Tap  on the My Glucose screen.
2. Select the tickbox next to the notes you would like to add. After you tick the box, you can add more specific information to your note.
 - Food notes: Enter meal type and grams or portion information
 - Insulin notes: Enter the number of units taken
 - Exercise notes: Enter intensity and duration
3. Tap **DONE** to save your note.

Notes you add are shown on your glucose graph and in your Logbook as symbols. You can review a note by tapping its symbol on your glucose graph or by going to the Logbook. See [Reviewing Your History](#) for more information about the Logbook. To edit a note from the glucose graph, tap the symbol and then tap on the information you would like to change. Tap **DONE** when you are finished.

Food

Insulin (Rapid or Long-acting)

Exercise

Food + insulin

Multiple/Custom notes – indicates different types of notes entered together or notes entered within a short period of time. A numbered badge next to the symbol indicates the number of notes.

Reviewing Your History

Reviewing and understanding your glucose history can be an important tool for improving your glucose control. The App stores about 90 days of information and has several ways to review your past glucose readings and notes. From the Main Menu, tap **Logbook** to view the Logbook or tap on one of the other history options under **Reports**.

IMPORTANT:

- Work with your health care professional to understand your glucose history.
 - Remember that FreeStyle LibreLink and Readers do not share data. For complete information on a device, be sure to scan your Sensor every 8 hours with that device; otherwise, your reports will not include all your data.
-

Logbook


The Logbook contains entries for each time you scanned your Sensor as well as notes you added. If you would like to view a different day, tap the  symbol or use the arrows. To add a note to a Logbook entry, tap on the entry and then tap . Select your note information and tap **DONE**.

To add a note that is independent of a Logbook entry, tap  on the main Logbook screen. Tap  if you want to add a note on a different date.


Other History Options

Daily Patterns: A graph showing the pattern and variability of your Sensor glucose readings over a typical day. The thick black line shows the median (midpoint) of your glucose readings. The light blue shading represents the 10th - 90th percentile range of your glucose readings. Dark blue shading represents the 25th - 75th percentile range.


Note: Daily Patterns needs at least 5 days of glucose data.

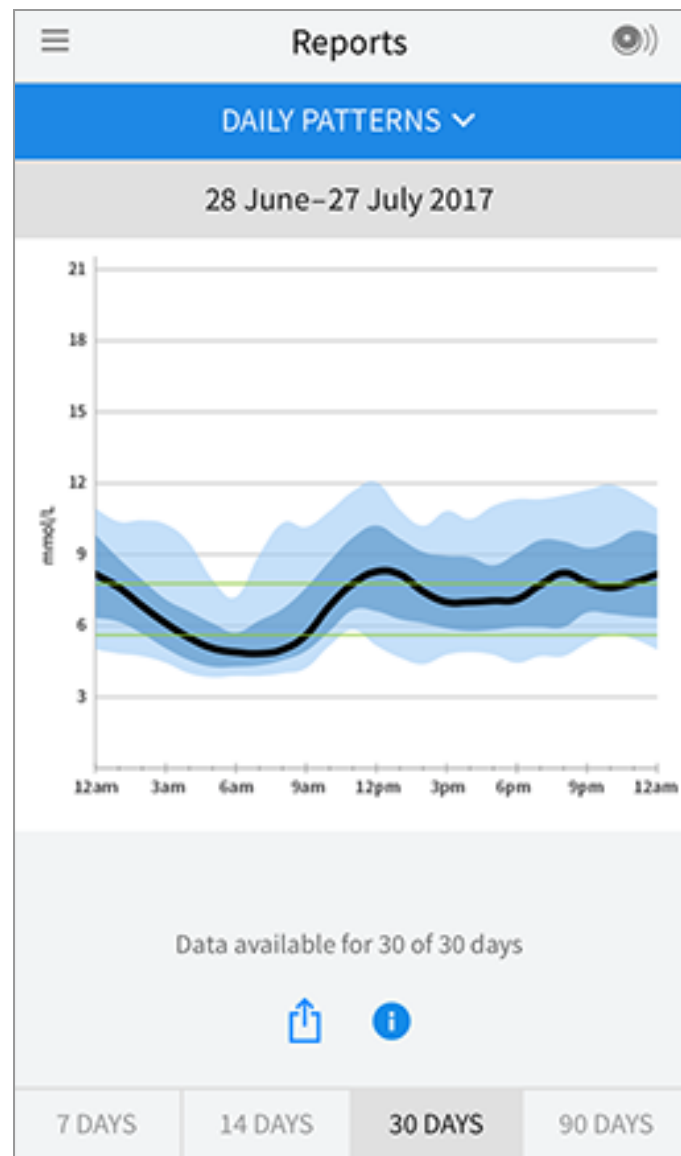
Time In Target: A graph showing the percentage of time your Sensor glucose readings were above, below, or within your Target Glucose Range.

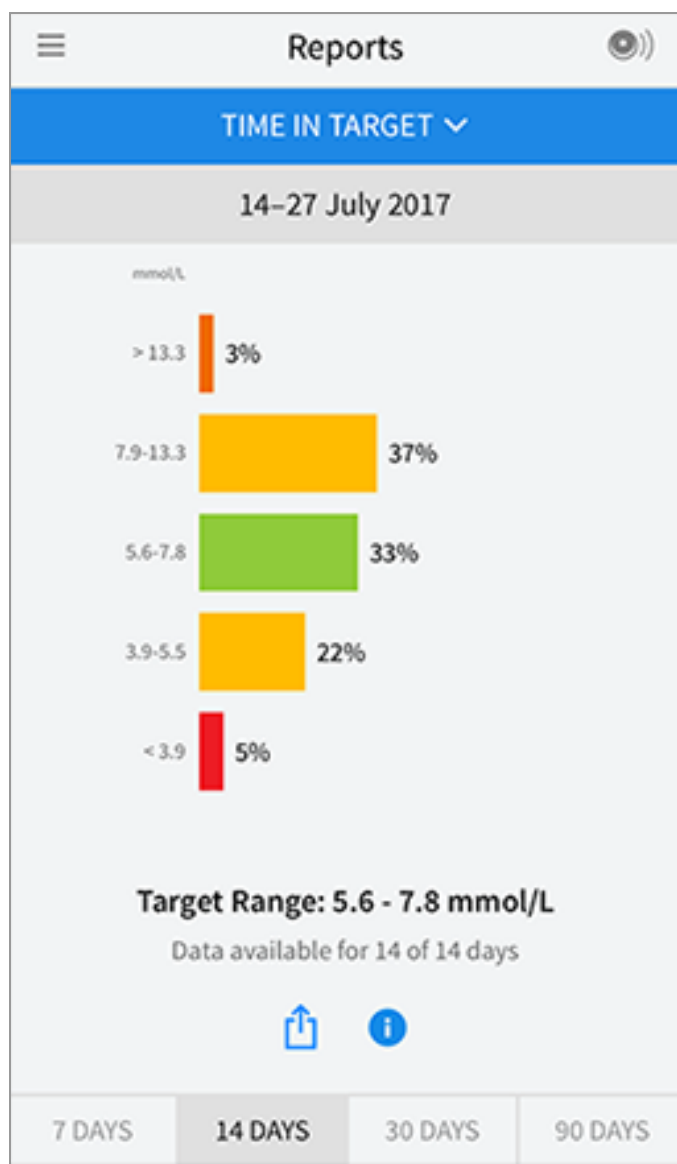Low Glucose Events: Information about the number of low glucose events measured by your Sensor. A low glucose event is recorded when your Sensor glucose reading is lower than 3.9 mmol/L for longer than 15 minutes. The total number of events is displayed below the graph. The bar graph displays the low glucose events in different periods of the day.

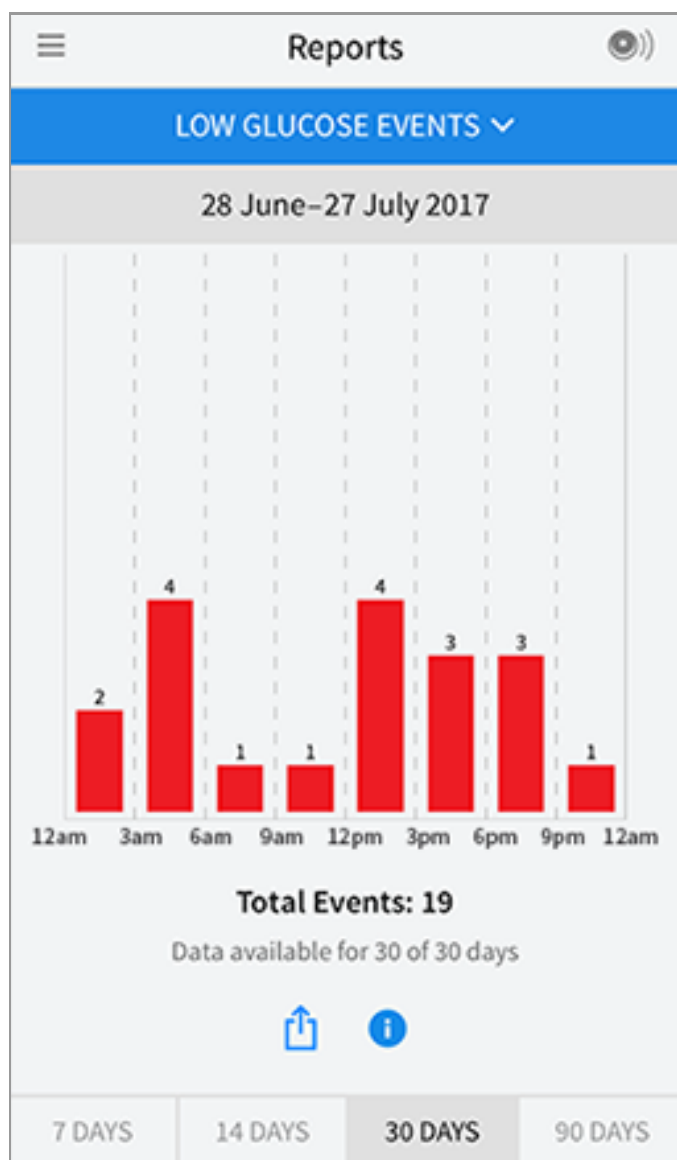
Average Glucose: Information about the average of your Sensor glucose readings. The overall average for the selected time period is displayed below the graph. The average is also shown for different periods of the day. Readings above or below your Target Glucose Range are yellow, orange, or red. Readings in range are green.

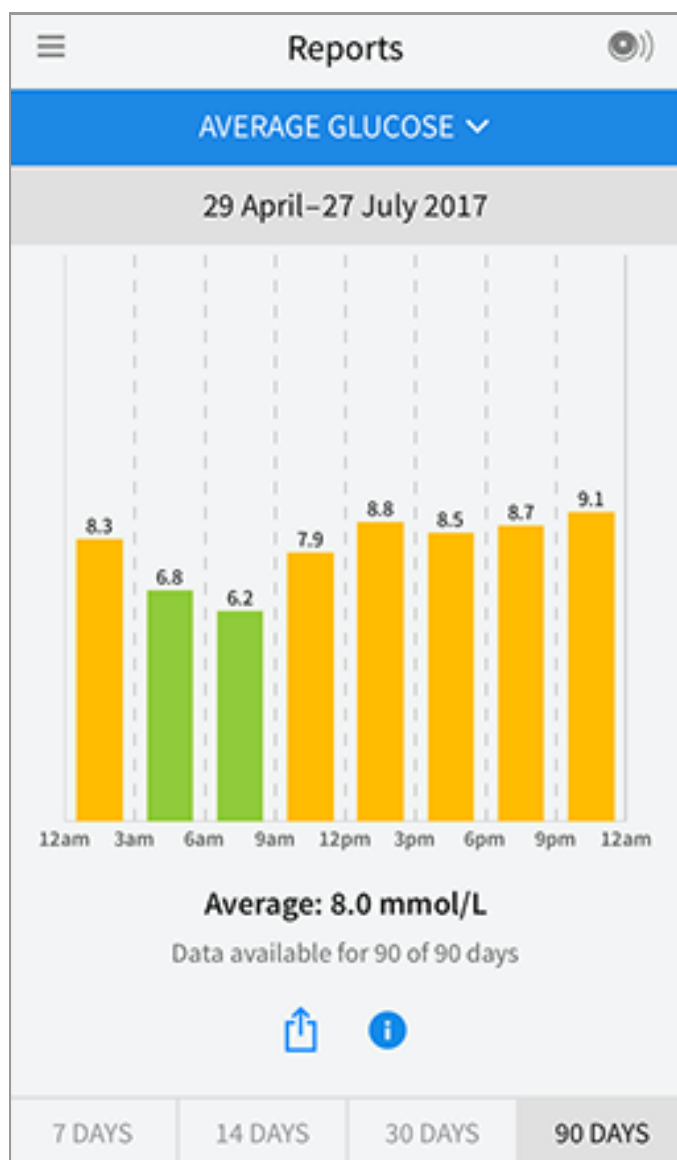Daily Graph: A graph of your Sensor glucose readings by day. The graph shows your Target Glucose Range and symbols for notes you have entered.

- The graph will scale to 27.8 mmol/L to accommodate glucose readings above 21 mmol/L.
- You might see gaps in the graph during times when you have not scanned at least once in 8 hours.
- The ⌚ symbol may appear indicating a time change. Gaps in the graph may result or glucose readings may be hidden.

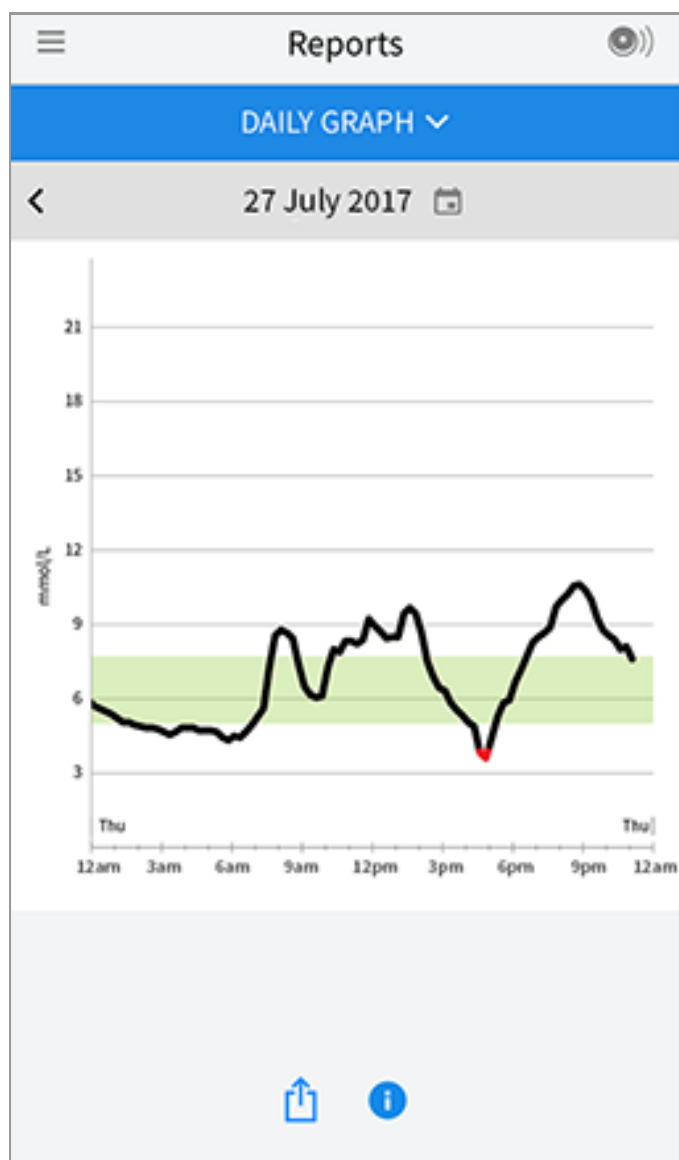
Estimated A1c: Your estimated A1c level (also called HbA1c) is based on available Sensor glucose data from the last 90 days. The more data available, the better your estimation will be. However, the estimated level may not match your A1c measured in a laboratory*. A1c can be used as an indicator of how well your glucose levels have been controlled and may be used to monitor your diabetes treatment regimen.


* The formula is based on the published reference, which compared average Sensor glucose and laboratory-measured A1c:

$$A1c_{\%} = (\text{Avg } SG_{\text{mmol/L}} + 2.59) / 1.59$$

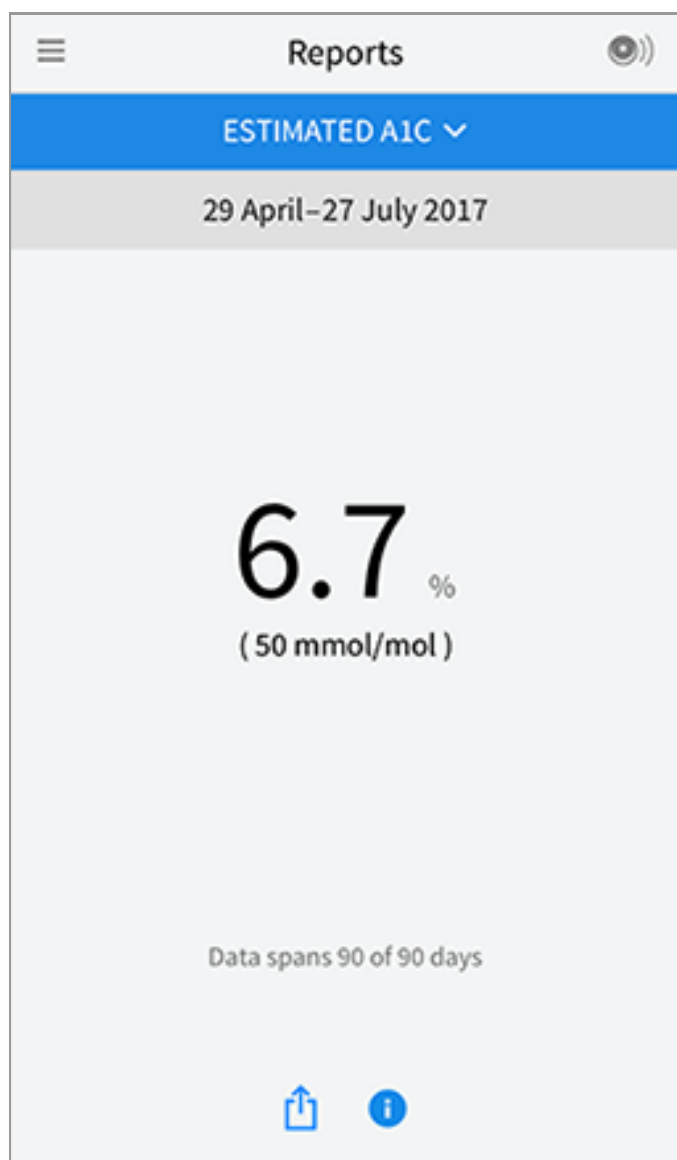Reference: Nathan DM, Kuenen J, Borg R, Zheng H, Schoenfeld D, Heine RJ for the A1c-Derived Average Glucose (ADAG). Study Group: Translating the hemoglobin A1c assay into estimated average glucose values. Diabetes Care 2008, 31:1473-8.

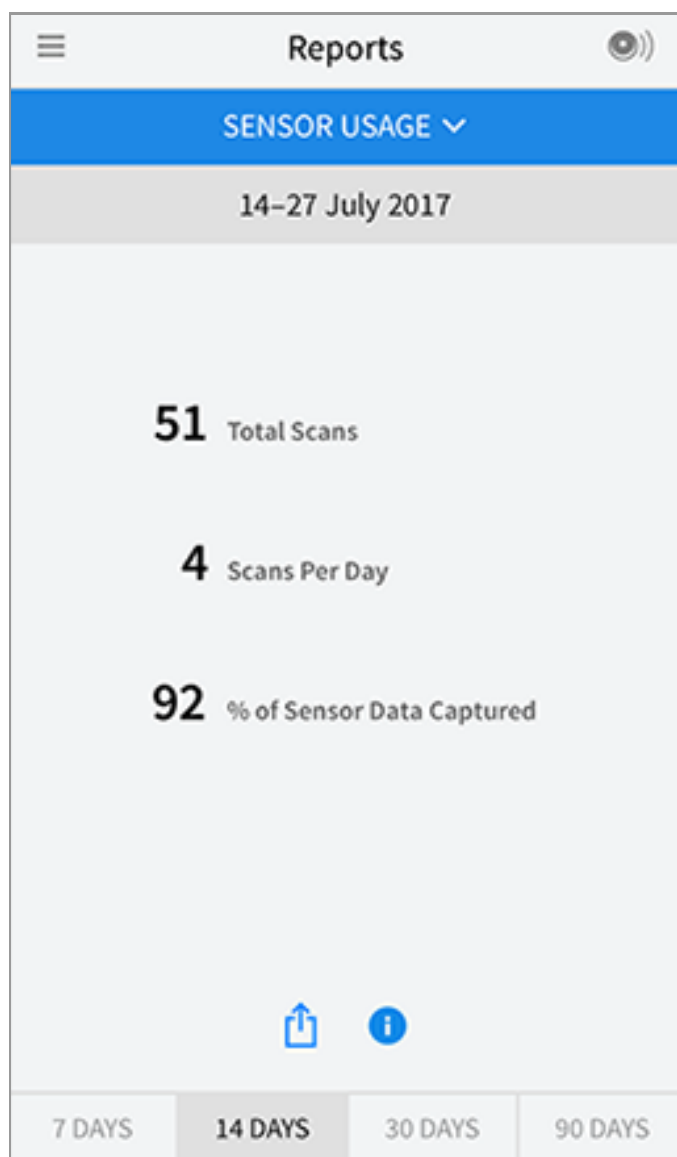Sensor Usage: Information about how often you scan your Sensor. This includes the total number of scans, an average of how many times you scanned your Sensor each day, and the percentage of possible Sensor data recorded from your scans.

Note:

- Tap the  symbol on any report to share a screenshot of the report.
- Tap the  symbol to view a description of the report.
- To view a different report, tap the drop down menu above the report, or go to the Main Menu.
- On all reports except the Daily Graph and Estimated A1c, you can select to show information about your last 7, 14, 30, or 90 days.

Removing Your Sensor

1. Pull up the edge of the adhesive that keeps your Sensor attached to your skin. Slowly peel away from your skin in one motion.

Note: Any remaining adhesive residue on the skin can be removed with warm soapy water or isopropyl alcohol.

2. Discard the used Sensor according to local regulations. See [Maintenance and Disposal](#). When you are ready to apply a new Sensor, follow the instructions in [Applying Your Sensor](#) and [Starting Your Sensor](#). If you removed your last Sensor before 14 days of use, you will be prompted to confirm that you would like to start a new Sensor when you first scan it.

Replacing Your Sensor

Your Sensor automatically stops working after 14 days of wear and must be replaced. You should also replace your Sensor if you notice any irritation or discomfort at the application site or if the App reports a problem with the Sensor currently in use. Taking action early can keep small problems from turning into larger ones.

CAUTION: If the glucose readings from the Sensor do NOT seem to match with how you feel, check to make sure that your Sensor has not come loose. If the Sensor tip has come out of your skin, or your Sensor is coming loose, remove the Sensor and apply a new one.

Setting Reminders

You can create single or repeating reminders to help you remember things like checking your glucose or taking insulin. There is one default reminder to help you remember to scan your Sensor. This Scan Sensor reminder can be changed or disabled but cannot be deleted.

Note: To receive reminders, make sure notifications for the App are enabled. If you want to receive a sound/vibration with your reminder, ensure that

sound/vibration on your smartphone is turned on, sound is set at a level you can hear, and your smartphone's Do Not Disturb feature (if available) is turned off. If Do Not Disturb is on, you will only see your reminder on the screen.

1. To add a new reminder, go to the Main Menu and tap **Reminders**. Tap **ADD REMINDER**.
2. Name your reminder.
3. Tap the time fields to set the time for the reminder.

Add Reminder

Reminder Name

Exercise

257

358

459 AM

500 PM

601

702

803

Repeating

☐ All

☐ Sunday

☒ Monday

☐ Tuesday

☒ Wednesday

☐ Thursday

☒ Friday

☐ Saturday

CANCEL

DONE

Note: If you would like the reminder to repeat, tap the slider to the right. You can also select which days you would like to receive the reminder.

4. Tap **DONE**. You will now see your reminder on the list along with the time you will receive it.

Note:

- To turn off a reminder, tap the slider to the left.
- To delete a reminder, swipe the reminder to the left and tap the  symbol.

The Scan Sensor reminder cannot be deleted.

- Your reminders will be received as notifications that you can swipe or tap to dismiss.

Settings and Other Options in the Main Menu

You can go to the Main Menu to change settings like your Target Glucose Range or LibreView password. You can also access the Share option, Help, and information about the App.

Settings

App Settings:

Unit of Measurement - View the glucose unit of measure used in the App.

Target Glucose Range - Set the target range you want to be displayed on the App glucose graph. It is also used to calculate your Time In Target. Tap **SAVE** when you are done.

Carbohydrate Units - Choose grams or portions for food notes that you enter. Tap **SAVE** when you are done.

Text to Speech - Turn on Text to Speech to have the glucose reading read aloud when you scan the Sensor. You will hear only your current glucose value and trend arrow direction along with whether you need to do a blood glucose test before making treatment decisions. Additional information, such as the glucose graph and any message, are available on your My Glucose screen. Always review your My Glucose screen to get complete information. Remember that Text to Speech inherits the volume settings on your smartphone. If your smartphone volume is turned off, you will not hear the glucose reading read aloud. Tap **SAVE** when you are done.

Account Settings:

Account Settings - View/change your LibreView account information.

Account Password - Change your LibreView account password.

Share

The Share option in the Main menu opens a web browser within the App. It lists different apps you can connect with to share your data. To connect your data with apps listed in the Share option, select them from the list of apps, and follow the onscreen instructions.

Help

View in-app tutorials, access this User's Manual, and review the App's legal information. You can also view the Event Log, which is a list of events recorded by the App. This may be used by Customer Service to help troubleshoot.

About

View App software version and other information.

Activities

Bathing, Showering, and Swimming: Your Sensor is water-resistant and can be worn while bathing, showering, or swimming. Do NOT take your Sensor deeper than 1 metre (3 feet) or immerse it longer than 30 minutes in water.

Sleeping: Your Sensor should not interfere with your sleep. It is recommended that you scan your Sensor before going to sleep and when you wake up because your Sensor holds 8 hours of data at a time. If you have reminders set to go off while you are sleeping, place your smartphone nearby.

Travelling by Air: Check with the airline prior to departure as rules and regulations may change without notice. Notify security personnel of the presence of your device when going through security systems. After you have put your smartphone in airplane mode, you can continue to get Sensor glucose readings. Do not scan your Sensor if restricted by flight regulations.

Note: Changing the time affects the graphs, statistics, and settings programmed by time of day since these are dependent on the time setting being accurate. The  symbol may appear on your glucose graph indicating a time change. Gaps in the graph may result or glucose readings may be hidden.

Maintenance and Disposal

Maintenance: The Sensor has no serviceable parts.

Disposal: The product contains electronic equipment, batteries, sharps and materials that may contact bodily fluids during use. Dispose of product in accordance with all applicable local regulations. Contact Customer Service for further information on the appropriate disposal of components.

Troubleshooting

This section lists problems that you may experience, the possible cause(s), and recommended actions. If there is an error, a message will appear on the screen with directions to resolve the error.

IMPORTANT: If you are having issues with the App, please keep in mind that uninstalling the App will cause you to lose all historical data and end the Sensor currently in use. Please call Customer Service if you have any questions.

Problems at the Sensor Application Site

Problem: **The Sensor is not sticking to your skin.**

What it may mean: The site is not free of dirt, oil, hair, or sweat.

What to do: 1. Remove the Sensor. 2. Consider shaving and/or cleaning the site with soap and water. 3. Follow the instructions in [Applying Your Sensor](#) and [Starting Your Sensor](#).

Problem: **Skin irritation at the Sensor application site.**

What it may mean: Seams or other constrictive clothing or accessories causing friction at the site **OR** you may be sensitive to the adhesive material.

What to do: Ensure that nothing rubs on the site. If the irritation is where the adhesive touches skin, contact your health care professional to identify the best solution.

Problems Starting Your Sensor or Receiving Sensor Readings

Display: **Sensor Starting Up**

What it may mean: Sensor is not ready to read glucose.

What to do: Wait until the 60 minute Sensor start-up period has completed.

Display: **Sensor Ended**

What it may mean: The Sensor life has ended.

What to do: Apply and start a new Sensor.

Display: New Sensor Found

What it may mean: You scanned a new Sensor before your previous Sensor ended.

What to do: Your smartphone can only be used with one Sensor at a time. If you start a new Sensor, you will no longer be able to scan your old Sensor. If you would like to begin using the new Sensor, select “Yes”.

Display: Sensor Error

What it may mean: The Sensor is unable to provide a glucose reading.

What to do: Scan again in 10 minutes.

Display: Glucose Reading is Unavailable

What it may mean: The Sensor is unable to provide a glucose reading.

What to do: Scan again in 10 minutes.

Display: Sensor Too Hot

What it may mean: Your Sensor is too hot to provide a glucose reading.

What to do: Move to a location where the temperature is appropriate and scan again in a few minutes.

Display: Sensor Too Cold

What it may mean: Your Sensor is too cold to provide a glucose reading.

What to do: Move to a location where the temperature is appropriate and scan again in a few minutes.

Display: Check Sensor

What it may mean: The Sensor tip may not be under your skin.

What to do: Try to start your Sensor again. If you see “Check Sensor” again on the screen, your Sensor was not applied properly. Apply and start a new Sensor.

Display: Replace Sensor

What it may mean: The App has detected a problem with your Sensor.

What to do: Apply and start a new Sensor.

Display: **Unexpected Application Error**

What it may mean: The App has detected an unexpected error.

What to do: Shut down the App completely and restart.

Display: **Incompatible Sensor**

What it may mean: The Sensor cannot be used with the App.

What to do: Call Customer Service.

Display: **Scan Error**

What it may mean: iPhone was unable to scan the Sensor.

What to do: Your scan was unsuccessful. Tap the scan button and scan again.

Customer Service

Customer Service is available to answer any questions you may have about FreeStyle LibreLink. Please go to www.FreeStyleLibre.com or refer to the product insert in your Sensor Kit for your Customer Service phone number. A printed copy of this User's Manual is available upon request.

Sensor Specifications

Sensor glucose assay method: Amperometric electrochemical sensor

Sensor glucose reading range: 2.2 to 27.8 mmol/L

Sensor size: 5mm height and 35mm diameter

Sensor weight: 5 grams

Sensor power source: One silver oxide battery

Sensor life: Up to 14 days

Sensor memory: 8 hours (glucose readings stored every 15 minutes)

Operating temperature: 10°C to 45°C

Sensor Applicator and Sensor Pack storage temperature: 4°C to 25°C

Operating and storage relative humidity: 10-90%, non-condensing

Sensor water resistance and ingress protection: IP27: Can withstand immersion into one metre (3 ft) of water for up to 30 minutes. Protected against insertion of objects >12mm diameter.

Operating and storage altitude: -381 metres (-1,250 ft) to 3,048 metres (10,000 ft)

Labelling Symbols and Definitions

 Consult instructions for use

 Temperature limit

 Manufacturer

 Batch code

 Type BF applied part

 Sensor code

 Do not re-use

 Use-by date

 Catalogue number

 Serial number

 Caution

 Sterilised by irradiation

Humidity limitation

Do not use if package is damaged

This product must not be disposed of via municipal waste collection. Contact the manufacturer for details.

Electromagnetic Compatibility

- The Sensor needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in this manual.
- Portable and mobile RF communications equipment can affect the Sensor.
- The use of accessories, transducers and cables other than those specified by Abbott Diabetes Care may result in increased EMISSIONS or decreased IMMUNITY of the Sensor.
- The Sensor should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the Sensor should be observed to verify normal operation in the configuration in which it will be used.
- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- Changes or modifications not approved by Abbott could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio

interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication. This device complies with Industry Canada's licence-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Guidance and manufacturer's declaration – electromagnetic emissions

The Sensor is intended for use in the electromagnetic environment specified below. The customer or the user of the Sensor should assure that it is used in such an environment.

Emissions test: RF emissions; CISPR 11

Compliance: Group 1

Electromagnetic environment – guidance: The Sensor uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.

Emissions test: RF emissions; CISPR 11

Compliance: Class B

Electromagnetic environment – guidance: The Sensor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.

Guidance and manufacturer's declaration – electromagnetic immunity

The Sensor is intended for use in the electromagnetic environment specified below. The customer or the user of the Sensor should assure that it is used in such an environment.

Immunity test: Electrostatic discharge (ESD); IEC 61000-4-2

IEC 60601 test level: ± 8 kV contact; ± 15 kV air

Compliance level: ± 8 kV contact; ± 15 kV air

Electromagnetic environment – guidance: Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.

Immunity test: Power frequency (50/60 Hz); magnetic field; IEC 61000-4-8

IEC 60601 test level: 30 A/m

Compliance level: 30 A/m

Electromagnetic environment – guidance: Power frequency magnetic fields should be at levels characteristic of a typical location in a typical domestic, commercial, or hospital environment.

Immunity test: Radiated RF; IEC 61000-4-3

IEC 60601 test level: 10 V/m; 80 MHz to 2.7 GHz

Compliance level: 10 V/m

Electromagnetic environment – guidance:

Recommended separation distance

$$d = 1.2 \sqrt{P}$$

80 MHz to 800 MHz

$$d = 2.3 \sqrt{P}$$

800 MHz to 2.5 GHz

P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).

Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b

Interference may occur in the vicinity of equipment marked with the following symbol:

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Sensor is used exceeds the applicable RF compliance level above, the Sensor should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Sensor.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the Sensor

The Sensor is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Sensor can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Sensor as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the

transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Performance Characteristics

Note: Please consult your health care team on how to use the information in this section. Refer to the Quick Reference Guide for a plain language summary of key performance measures.

Clinical Study Overview

Performance of the FreeStyle Libre Flash Glucose Monitoring System (the System) was evaluated in a clinical study.

The study was conducted at 4 centres with a total of 95 subjects with diabetes (84.2% Type 1, 15.8% Type 2). All subjects were aged eighteen and older. Subjects in the study required insulin to manage their diabetes. Each subject wore up to two FreeStyle Libre Sensors on the back of the upper arm. During the study, subjects tested their blood glucose using fingerstick capillary samples at least eight times during each day of the study. Subjects used the blood glucose meter built into the FreeStyle Libre Reader. Additionally, subjects had their venous blood glucose analysed up to 112 times over four separate visits to the clinical centre. Venous blood was analysed using the Yellow Springs Instrument Life Sciences 2300 STAT Plus™ Glucose & Lactate Analyzer (YSI). YSI is a laboratory glucose and lactate analyser of whole blood and plasma and is a widely recognised standard in laboratory analysis of blood glucose. Glucose readings obtained from the System were compared to glucose readings obtained from the YSI to evaluate the performance of the System. Three lots of Sensors were evaluated in the study.

Agreement with YSI Levels

Agreement between FreeStyle Libre Glucose Measurement (GM) and venous blood was characterised by using paired GM and Yellow Springs Instrument measurements (YSI). The accuracy of GM versus YSI reference was assessed by calculating the percentage of System readings that were within 15%, 20%, 30% and 40% for reference values 4.44 mmol/L and above, and 0.83 mmol/L, 1.11 mmol/L, 1.67 mmol/L and 2.22 mmol/L for values below 4.44 mmol/L when glucose levels are assigned using the YSI values. Overall, 90.3% of results were within ± 1.11 mmol/L / 20% of YSI reference.

Agreement with GM Glucose Levels

Agreement between GM and venous blood was characterised by using paired GM and Yellow Springs Instrument measurements (YSI). The accuracy of GM versus YSI reference was assessed by calculating the percentage of System readings that were within 15%, 20%, 30% and 40% for reference values 4.44 mmol/L and above, and 0.83 mmol/L, 1.11 mmol/L, 1.67 mmol/L and 2.22 mmol/L for values below 4.44 mmol/L. The results are presented in **Table 1** for YSI reference. Overall, 90.5% of results were within ± 1.11 mmol/L / 20% of YSI reference.

Table 1: Number and Percent of Results within YSI Reference

GM Glucose Level (mmol/L)	Number of GM-Reference Pairs	Within ±15% / ±0.83 mmol/L	Within ±20% / ±1.11 mmol/L	Within ±30% / ±1.67 mmol/L	Within ±40% / ±2.22 mmol/L	Outside ±40% / ±2.22 mmol/L
Overall	9766	82.6	90.5	96.9	99.0	1.0
2.2-2.8	43	48.8	58.1	74.4	86.0	14.0
2.8-4.4	767	55.1	67.3	85.1	93.9	6.1
4.5-10.0	6120	82.2	90.7	97.4	99.4	0.6
10.0-16.7	2517	91.6	96.8	99.4	99.8	0.2
16.7-22.2	291	89.0	95.2	99.7	99.7	0.3
22.3-27.8	28	100.0	100.0	100.0	100.0	0.0

Agreement on Day 1 against YSI Reference

The accuracy of GM versus YSI reference on the first day of Sensor wear was assessed by calculating the percentage of System readings that were within 15%, 20%, 30% and 40% for reference values 4.44 mmol/L and above, and within 0.83 mmol/L, 1.11 mmol/L, 1.67 mmol/L and 2.22 mmol/L for values below 4.44 mmol/L by hourly intervals. The results are presented in **Table 2**.

Table 2: Number and Percent of Results within YSI Reference

Time Interval (hours)	Number of GM-Reference Pairs	Within ±15% / ±0.83 mmol/L	Within ±20% / ±1.11 mmol/L	Within ±30% / ±1.67 mmol/L	Within ±40% / ±2.22 mmol/L	Outside ±40% / ±2.22 mmol/L
(0-2] [*]	182	76.9	87.4	95.1	98.9	1.1
(2-4]	385	70.6	82.9	94.3	98.4	1.6
(4-6]	381	76.9	85.3	91.3	96.1	3.9
(6-8]	381	75.9	85.3	94.8	98.2	1.8
(8-16] [†]	107	73.8	83.2	97.2	100.0	0.0
(16-18]	303	80.2	89.4	98.3	100.0	0.0
(18-20]	344	80.8	87.5	95.3	99.4	0.6
(20-22]	336	81.5	89.6	97.9	98.8	1.2
(22-24]	155	86.5	93.5	97.4	100.0	0.0

* (0-2] interval includes the 60 minute start-up time.

† Time interval is not divided to smaller segments due to small number of data pairs.

Overall Accuracy against YSI reference

Accuracy was measured by comparing the absolute relative difference between the System and reference YSI glucose values. The absolute relative difference measures the level of disagreement between the System and the reference value, but does not tell you whether the System glucose value was, on average, higher or lower than the reference glucose value. The Mean Absolute Relative Difference gives an indication of the average percent disagreement between the GM and the reference. **Table 3** shows the overall absolute difference measure. Overall, the Mean Absolute Relative Difference was 9.5% for the comparison with YSI reference. The Median Absolute Relative Difference shows that half of the time the System was within 7.3% of the YSI reference.

Table 3: Difference Measures with YSI Reference

Number of GM-Reference Pairs	Median Absolute Relative Difference (%)	Mean Absolute Relative Difference (%)
9766	7.3	9.5

Agreement with BG Levels

Agreement between the System and capillary blood glucose values (BG) as measured by the Reader’s built-in meter was characterised by using paired System GM and BG value. The accuracy of GM versus BG value was assessed by calculating the percentage of System readings that were within 15%, 20%, 30% and 40% for BG values 4.44 mmol/L and above, and within 0.83 mmol/L, 1.11 mmol/L, 1.67 mmol/L and

2.22 mmol/L for values below 4.44 mmol/L. The results are presented in **Table 4** for BG values. Overall, 86.9% of results were within ± 1.11 mmol/L / 20% of BG values.

Table 4: Number and Percent of Results within BG Values*

GM Glucose Level (mmol/L)	Number of GM-Reference Pairs	Within $\pm 15\%$ / ± 0.83 mmol/L	Within $\pm 20\%$ / ± 1.11 mmol/L	Within $\pm 30\%$ / ± 1.67 mmol/L	Within $\pm 40\%$ / ± 2.22 mmol/L	Outside $\pm 40\%$ / ± 2.22 mmol/L
Overall	9250	75.9	86.9	95.9	98.5	1.5
2.2-2.8	127	52.0	59.8	74.8	89.0	11.0
2.8-4.4	807	59.9	70.8	86.5	93.8	6.2
4.5-10.0	4841	73.7	86.0	96.3	99.0	1.0
10.0-16.7	2798	83.7	92.7	98.4	99.4	0.6
16.7-22.2	563	82.4	93.6	98.0	99.3	0.7
22.3-27.8	114	80.7	93.0	97.4	98.2	1.8

* Comparison to BG was performed using the Reader's built-in blood glucose meter. Different performance may be expected when compared to other models of blood glucose meters.

Overall Accuracy against BG values

Accuracy was measured by comparing the absolute relative difference between the System and BG values. The absolute relative difference measures the level of disagreement between the System and the BG value, but does not tell you whether the System glucose value was, on average, higher or lower than the BG glucose value. The Mean Absolute Relative Difference gives an indication of the average percent disagreement between the GM and the BG value. **Table 5** shows the overall absolute difference measure. Overall, the Mean Absolute Relative Difference was 11.4% for the comparison with BG value. The Median Absolute Relative Difference shows that half of the time the System was within 9.1% of the BG value.

Table 5: Difference Measures with BG Value*

Number of GM-Reference Pairs	Median Absolute Relative Difference (%)	Mean Absolute Relative Difference (%)
9250	9.1	11.4

* Comparison to BG was performed using the Reader's built-in blood glucose meter. Different performance may be expected when compared to other models of blood glucose meters.

Concurrence of System and Reference (GM vs. YSI)

The percentage of concurring glucose values (GM vs. YSI) in each glucose reference range is presented for each GM range in **Table 6**. For example, when the System glucose results were within the 4.5 to 6.7 mmol/L range, actual blood glucose values were less than 2.2 mmol/L 0.0% of the time, between 2.2 and 3.3 mmol/L 0.1% of the time, between 3.4 and 4.4 mmol/L 2.7% of the time, between 4.5 and 6.7 mmol/L 67.9% of the time, between 6.7 and 8.9 mmol/L 28.0% of the time, between 8.9 and 11.1 mmol/L 1.2% of the time, between 11.2 and 13.9 mmol/L 0.1% of the time and above 13.9 mmol/L 0.0% of the time.

Table 6: Concurrence Analysis by Glucose Level

GM (mmol/L)	YSI Glucose Level (mmol/L)												N
	<2.2*	2.2-3.3	3.4-4.4	4.5-6.7	6.7-8.9	8.9-11.1	11.2-13.9	13.9-16.7	16.7-19.4	19.5-22.2	22.3-27.8	>27.8*	
<2.2	0.0	42.1	36.8	21.1	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	19
2.2-3.3	2.6	28.8	45.5	23.1	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	156
3.4-4.4	0.2	4.9	32.7	60.2	2.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	654
4.5-6.7	0.0	0.1	2.7	67.9	28.0	1.2	0.1	0.0	0.0	0.0	0.0	0.0	2497
6.7-8.9	0.0	0.0	0.0	6.9	71.4	20.5	0.9	0.2	0.0	0.0	0.0	0.0	2635
8.9-11.1	0.0	0.0	0.0	0.2	11.0	68.8	18.5	1.1	0.3	0.0	0.0	0.0	1761
11.2-13.9	0.0	0.0	0.0	0.0	0.3	15.3	72.5	11.6	0.3	0.0	0.0	0.0	1190
13.9-16.7	0.0	0.0	0.0	0.0	0.0	0.5	23.5	65.3	10.6	0.0	0.0	0.0	554
16.7-19.4	0.0	0.0	0.0	0.0	0.0	0.0	0.0	33.3	56.6	7.3	0.0	0.0	219
19.5-22.2	0.0	0.0	0.0	0.0	0.0	0.0	0.0	1.4	20.8	66.7	11.1	0.0	72
22.3-27.8	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	21.4	78.6	0.0	28
>27.8	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	92.3	7.7	13

* Levels out of system dynamic range.

Concurrence Analysis by Glucose Trend Arrow

Trend Arrow accuracy of the System, as assessed by concurrence analysis, is presented in **Table 7**. For example, when the trend arrow indicated that glucose was changing slowly (-0.06 to 0.06 mmol/L/min (→)), actual glucose levels in the body were falling quickly (↓) 0.2% of the time, falling (↘) 2.6% of the time, changing slowly (→) 86.6% of the time, rising (↗) 4.2% of the time, and were rising quickly (↑) 0.6% of the time.

Table 7: Concurrence Analysis by Glucose Trend Arrow

GM (mmol/L /min)	YSI (mmol/L/min)							N
	<-0.1	[-0.1, -0.06]	[-0.06, 0]	[0, 0.06]	[0.06, 0.1]	>0.1	NA*	
<-0.1 (↓)	39.6	42.4	9.0	0.7	0.0	0.0	8.3	144
-0.1 to -0.06 (↘)	4.7	27.2	55.5	5.9	0.3	0.0	6.4	699
-0.06 to 0.06 (→)	0.2	2.6	50.2	36.4	4.2	0.6	5.8	7809
0.06 to 0.1 (↗)	0.0	0.6	5.7	43.4	36.1	9.8	4.4	723
>0.1 (↑)	0.0	0.0	2.0	12.9	36.9	43.0	5.2	249
NA†	0.6	8.0	28.2	37.4	12.6	4.0	9.2	174

* Glucose rate of change not available due to the time difference between glucose readings exceeding 30 minutes.

† Glucose Trend Arrow not available.

Agreement with ‘LO’ and ‘HI’ GM Reading against YSI Reference

The System reports glucose concentrations between 2.2 and 27.8 mmol/L. When the System determines that glucose level is below 2.2 mmol/L, it will report as ‘LO’. When the System determines that glucose level is above 27.8 mmol/L, it will report as ‘HI’. **Table 8** displays the concurrence between the GM and YSI reference glucose when GM reads ‘LO’. For example, when GM reading was ‘LO’, YSI glucose values were less than 2.2 mmol/L 0.0% of the time, equal to or above 2.2 mmol/L 100.0% of the time, above 2.8 mmol/L 84.2% of the time, above 3.3 mmol/L 57.9% of the time, above 3.9 mmol/L 36.8% of the time, and above 4.4 mmol/L 21.1% of the time.

Table 8: Concurrence Analysis with ‘LO’ GM Reading

% of GM points in YSI range	YSI (mmol/L)						N
	<2.2	≥2.2	>2.8	>3.3	>3.9	>4.4	
	0.0	100.0	84.2	57.9	36.8	21.1	
Number of GM points in YSI range	0	19	16	11	7	4	19

Table 9 displays the concurrence between the GM and YSI reference glucose when GM reads 'HI'. For example, when GM reading was 'HI', YSI glucose values were less than or equal to 11.1 mmol/L 0.0% of the time, above 11.1 mmol/L 100.0% of the time, above 16.7 mmol/L 100.0% of the time, above 22.2 mmol/L 100.0% of the time, and above 27.8 mmol/L 7.7% of the time.

Table 9: Concurrence Analysis with 'HI' GM Reading

% of GM points in YSI range	YSI (mmol/L)					N
	≤11.1	>11.1	>16.7	>22.2	>27.8	
	0.0	100.0	100.0	100.0	7.7	
Number of GM points in YSI range	0	13	13	13	1	13

Accuracy by Day of Wear

The Sensor can be worn for up to 14 days. To show Sensor performance over time, the absolute relative difference between the System and reference YSI glucose values over the wear duration is presented in **Table 10**.

Table 10: Difference Measures by Day: (YSI Reference)

Day	Number of GM-Reference Pairs	Median Absolute Relative Difference (%)	Mean Absolute Relative Difference (%)
1	2578	8.3	10.9
6	2544	7.1	8.8
11	2419	7.4	9.3
14	2225	6.6	9.0

The accuracy of GM versus YSI reference was assessed by calculating the percentage of System readings that were within 15%, 20%, 30% and 40% for reference values of 4.44 mmol/L and above, and 0.83 mmol/L, 1.11 mmol/L, 1.67 mmol/L and 2.22 mmol/L for values below 4.44 mmol/L. The results are presented in **Table 11**.

Table 11: Number and Percent of Results within YSI Reference

Day	Number of GM-Reference Pairs	Within ±15% / ±0.83 mmol/L	Within ±20% / ±1.11 mmol/L	Within ±30% / ±1.67 mmol/L	Within ±40% / ±2.22 mmol/L	Outside ±40% / ±2.22 mmol/L
1	2578	77.8	86.8	95.4	98.6	1.4
6	2544	86.0	92.9	98.2	99.4	0.6
11	2419	82.3	90.2	98.3	99.9	0.1
14	2225	84.0	91.4	96.3	99.1	0.9

System Glucose Availability

The System is designed to produce a glucose reading after each user-initiated scan that is performed after the start-up time, through the wear period. **Table 12** shows the number of available glucose readings reported by all Sensors and the expected number based on the total number of scan attempts. Results are shown for Sensors which produced at least one GM reading during the clinical study over the total wear period. The percentage of available GM readings is presented in comparison to the number of expected GM readings. Overall, 99.7% of GM readings (18606 GM readings out of an expected 18666) were available.

Table 12: GM Availability

No. GM	No. Scan	%
18606	18666	99.7

Detection of Hypoglycaemic and Hyperglycaemic Events

Table 13 shows the accuracy of the System's Glucose Messages in informing the user of low or high

glucose events.

Percentages are displayed for three different parameters:

- Detection Rate – amount of time the System displays a Glucose Message correctly.
- Missed Detection Rate – amount of time the System did not display a Glucose Message when it should have.
- False Notification Rate – amount of time the System displays a Glucose Message when it shouldn't have.

Low and High Glucose Message percentages are based on comparison between Sensor glucose and true blood glucose values within 15 minutes before or after the true blood glucose value.

Glucose Going Low and Glucose Going High Message percentages are based on comparison between Sensor glucose and true blood glucose values up to 15 minutes before the true blood glucose value.

For example, the System was able to detect 77.6% of actual low glucose events (detection rate), but 54.7% of the time a Low Glucose message was displayed in error (false notification rate) and 22.4% of the time a Low Glucose message was not displayed when it should have been (missed detection rate).

Table 13: Detection of Hypoglycaemic and Hyperglycaemic Events

Type of Notification	Notification Status	15 Minute Interval
Notification of Hypoglycaemic Events (Low Glucose message)	Detection Rate (%)	77.6
	Missed Detection Rate (%)	22.4
	False Notification Rate (%)	54.7
Notification of Hyperglycaemic Events (High Glucose message)	Detection Rate (%)	84.9
	Missed Detection Rate (%)	15.1
	False Notification Rate (%)	10.6
Impending Notification of Hypoglycaemic Events (Glucose Going Low message)	Detection Rate (%)	83.2
	Missed Detection Rate (%)	16.8
	False Notification Rate (%)	70.0
Impending Notification of Hyperglycaemic Events (Glucose Going High message)	Detection Rate (%)	83.5
	Missed Detection Rate (%)	16.5
	False Notification Rate (%)	20.5

Precision

Precision of the System was evaluated by comparing the results from two separate Sensors worn on the same subject at the same time. **Table 14** provides data from two separate Sensors worn on 95 subjects at the same time. Overall 15300 real-time pairs of GM measurements were obtained, with a mean CV of 5.7%.

Table 14: Overall between Sensor Precision

Mean Glucose (mmol/L)	Median CV	Mean CV	Number of Subjects	Number of Paired Readings
9.0	4.2	5.7	95	15300

Sensor Wear Duration

The Sensor can be worn for up to 14 days. To estimate how long a Sensor will work over the wear duration, 95 Sensors were evaluated to determine how many days of readings each Sensor provided. Of the 95 Sensors, 82 (86.3%) lasted until the final day of use. 93 Sensors (97.9%) lasted more than 5 days.

There were 13 (13.7%) Sensors that failed early, of which 2 (2.1%) failed on or before the fifth day of wear.

Adverse Events

No device related serious adverse events occurred during the study. Mild skin irritations, such as erythema, bleeding, bruising, infection and papule were reported around the insertion site and adhesive area by a moderate frequency of subjects (8 out of 95 or 8.4%).

Customer Service: www.FreeStyleLibre.com


Patent: <https://www.abbott.com/patents>

FreeStyle, Libre, and related brand marks are trademarks of Abbott Diabetes Care Inc. in various jurisdictions. Other trademarks are the property of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries.

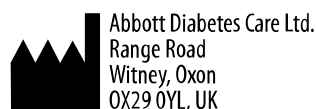App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

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©2019 Abbott
Apple et le logo Apple sont des marques commerciales d'Apple Inc., déposées aux États-Unis et dans d'autres pays.
App Store est une marque de service d'Apple Inc., déposée aux États-Unis et dans d'autres pays.