



Quick Start Guide for FreeStyle LibreLink app

A FreeStyle Libre product




Setup Overview

Refer to your User’s Manual for full system instructions and information.



iPhone
Download **FreeStyle LibreLink** from App Store
Android phone
Download **FreeStyle LibreLink** from Google Play Store

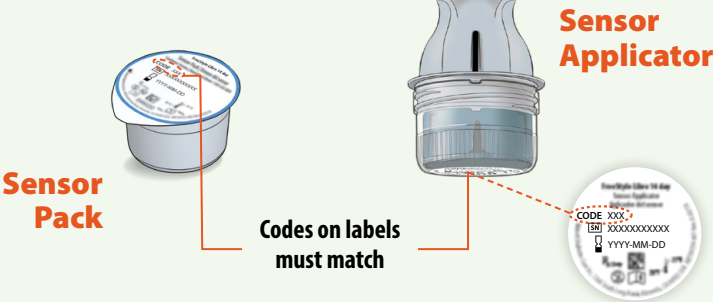


- 1 Assemble and Apply Sensor to the back of your upper arm
- 2 Set up App
- 3 Use the App to scan your Sensor

After the start-up period, you can check your glucose

1 Assemble and Apply Sensor to the back of your upper arm

Get Ready

STEP 1

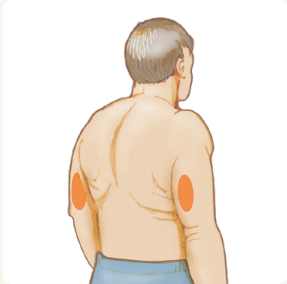


Select site on back of upper arm. Do not use other sites as these are not approved and may result in inaccurate glucose readings.
Note: Avoid scars, moles, stretch marks, lumps, and insulin injection sites. To prevent skin irritation, rotate sites between applications.

STEP 2



Clean site with alcohol wipe. Allow site to dry before proceeding.

STEP 3

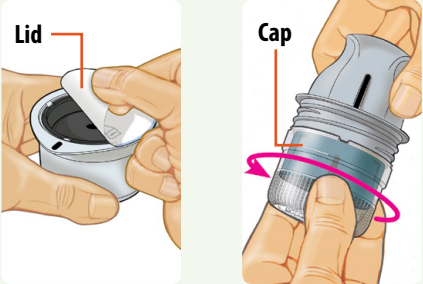


Peel lid completely off Sensor Pack. Unscrew cap from Sensor Applicator.
CAUTION: Sensor codes must match on Sensor Pack and Sensor Applicator or glucose readings will be incorrect.



Prepare Sensor Applicator

STEP 4

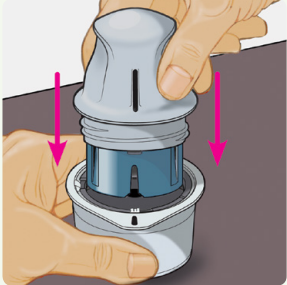

Line up dark mark on Sensor Applicator with dark mark on Sensor Pack. On a hard surface, press down firmly on Sensor Applicator until it comes to a stop.

STEP 5



Lift Sensor Applicator out of Sensor Pack.


STEP 6



Sensor Applicator is ready to apply Sensor.
CAUTION: Sensor Applicator now contains a needle. Do not touch inside Sensor Applicator or put it back into Sensor Pack.


Apply Sensor

STEP 7

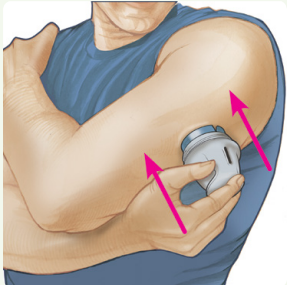


Place Sensor Applicator over site and push down firmly to apply Sensor.
CAUTION: Do not push down on Sensor Applicator until placed over prepared site to prevent unintended results or injury.

STEP 8

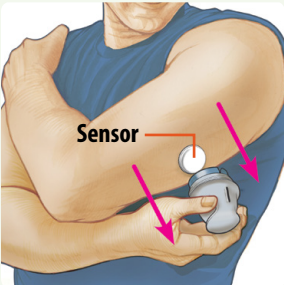

Gently pull Sensor Applicator away from your body.

STEP 9

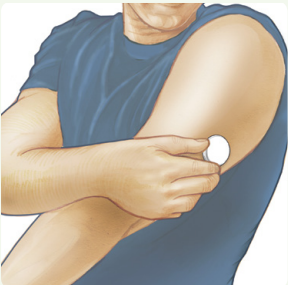


Make sure Sensor is secure. Discard used Sensor Applicator and Sensor Pack according to local regulations.



2 Set up App

STEP 1

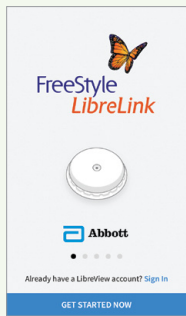

Check your smartphone is connected to a WiFi or cellular network. Tap **Sign In** or **GET STARTED NOW** .

Review the legal information and follow instructions to set up a LibreView account or login to your existing account.

STEP 2



Enter your personal settings. Tap **NEXT** after each screen.

STEP 3



Review important information about glucose readings, when to do a blood glucose test, and where to apply a Sensor. Tap **NEXT** to proceed until you see the Scan New Sensor screen.

3 Use the App to scan your Sensor

How to Scan

iPhone: Tap scan button (📶) . You can either tap the blue box on the Home Screen or (📶) at the top right. Hold top of iPhone near Sensor until you get a tone and/or vibration.

Android phone: Hold back of Android phone near Sensor until you get two tones and/or vibrations.


STEP 1



Scan your Sensor to start it.


Note: If you want to use both the Reader and the App to check your glucose, you must start the Sensor with the Reader first.

STEP 2

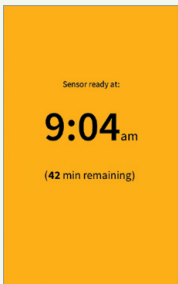

The Sensor can be used to check your glucose after the start-up period. While the Sensor is starting up, you can navigate away from the App.


STEP 3

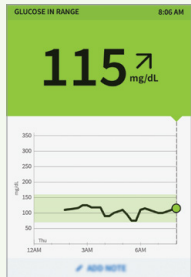

After the start-up period, scan your Sensor to get your glucose reading.

Glucose Graph


Graph of your current and stored glucose readings

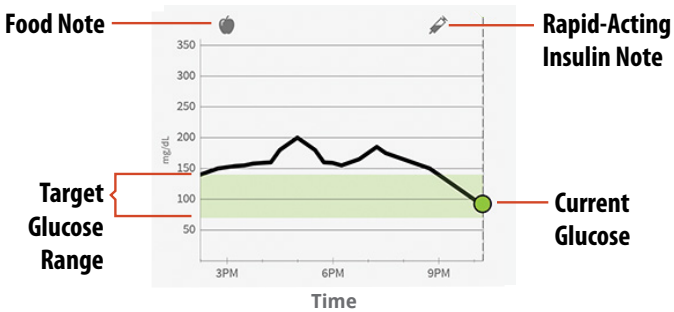

Message
Tap for more information

Check Blood Glucose
When you see this symbol, do a blood glucose test before making treatment decisions

Current Glucose
Glucose from your latest scan

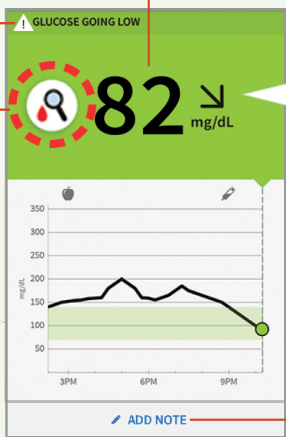

Add Note
Tap to add notes to glucose reading


Glucose Trend Arrow

Direction your glucose is going


Arrow	What it means
↑	Glucose rising quickly
↗	Glucose rising
→	Glucose changing slowly
↘	Glucose falling
↓	Glucose falling quickly

WARNING: When you see the 📶 symbol, you must check your blood glucose with a blood glucose meter before making any treatment decisions. Sensor readings may not accurately reflect blood glucose levels.


Reports of your glucose history

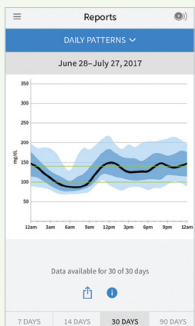


Daily Patterns

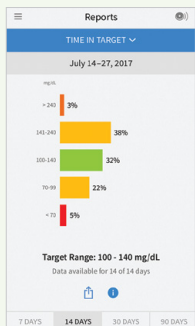


Time In Target

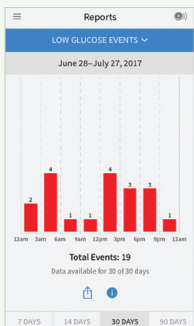


Low Glucose Events

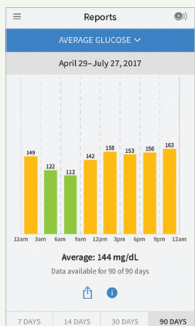

Average Glucose

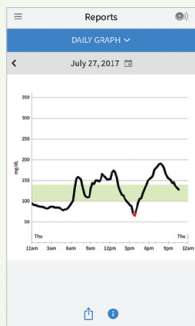

Daily Graph

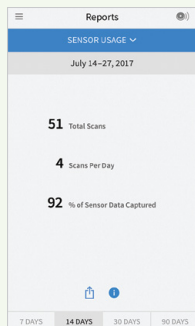

Sensor Usage

- Tap at the top left of the screen to access your reports.
- Tap (Android phone) or (iPhone) to share a screenshot of the report.
- Tap to view a description of the report.

Example screens displayed for iPhone

If you have followed the instructions described in the User's Manual and are still having difficulty setting up your System or if you are experiencing symptoms that are not consistent with your glucose monitoring, contact your health care professional.

FreeStyle, Libre, and related brand marks are trademarks of Abbott Diabetes Care Inc. in various jurisdictions.
©2019 Abbott ART39766-001 Rev. A 04/19

Consult instructions for use

Manufacturer
Abbott Diabetes Care Inc.
1360 South Loop Road
Alameda, CA 94502 USA

